



ZEBRONICS
ALWAYS AHEAD

Since 1997



Health Care Products

Zebronics has an extensive range of healthcare products that helps in effective health monitoring and sanitization. Our wide product range includes Contactless Infrared Thermometer, Disinfectant Fogging Machine, Portable Atomizer, Multi-function UV Sterilizer Box, Fingertip Pulse Oximeter, and PPE kits that help in keeping your homes, offices, schools, malls, and more such places sanitized and safe.

Page of Content

03	Reusable Outdoor Protection Mask with Nose Pin Z95	05	PPE Kit - Zeb-Isolation kit
07	PPE Kit - Zeb-Protection kit	09	Disinfectant Fogging Machine ZEB-FM500G
11	Disinfectant Fogging Machine	13	Contactless Infrared Thermometer AD801
15	Contactless IR Thermometer ZEB-DT100	17	Multifunction UV Sterilizer Box ZEB-UVS500
19	Portable UV Sterilizer ZEB-UVS200L	21	Fingertip Pulse Oximeter ZEB-FPO500
23	Portable Multi-purpose Atomizer ZEB-AM500P		

Reusable Outdoor Protection Mask with Nose Pin.

Z95



Reusable
(upto 30 gentle
hand washes)



Super
Breathable



Splash
Resistance



MRP


₹199/-



Anti-Dust



Anti-Bacterial

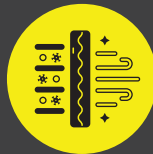

Anti-Pollution



6 layer
Filtration

Reusable Outdoor Protection Mask with Nose Pin.

Z95

Features

- Comfortable to wear and breathe
- Wider size and maximum face coverage
- Innovative 2 panel design for universal fit
- Adjustable nose clip for enhanced protection
- Comfortable and soft, elastic ear loops
- Minimal eyewear fogging
- Convenient flat-fold design for easy storage



Instructions for Use



Precautions

- Wash hands properly before wearing and removing the mask
- Wash the mask after everyday use, discard after 30 washes
- Store in dry place & away from sunlight
- Do not share the mask with anyone, for individual use only
- Change the mask if there is difficulty in breathing


Warning

- Improper usage may cause health risks
- Instructions for use to be followed strictly
- Mask should be discarded properly

Folding Instruction



Washing Instruction

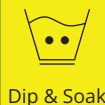

Dip & Soak



Do Not Bleach



Do Not Tumble Dry



Do Not Iron



Do Not Dry Clean



Do Not Wring



Drip Dry



ZEB- ISOLATION KIT

ISOLATION SUIT x 1



- 100gsm poly laminated fabric.
- Universal fit for all size, comfortable to wear.
- Made with finest needles of 12mm to avoid needle holes which prevents Bacteria/Virus/ Blood to penetrate.

MRP
₹1,450/-





SHOE COVER x 1 PAIR

Maximum Height for better protection.



EYE GOGGLES X1

Complete coverage, end to end protection.



FACE SHIELD X 1

Anti fog treated, Foam for comfortable headrest, Elastic fastener to fit every head size.



NITRILE GLOVES x 1 PAIR

High quality powder free nitrile gloves.



3 PLY MASK X 2

3 Layers with anti bacterial filter and nose pin for better comfort.



BIODEGRADABLE DISPOSABLE BAG X 1

All the items above after use can be collected in this bio hazard bag and discarded.



**COVID19/INMAS/
140520-17/RATL-01**

ZEB- PROTECTION KIT



PROTECTION SUIT x 1



- 🛡️ 90gsm poly laminated fabric.
- 🛡️ Universal fit for all size, comfortable to wear.
- 🛡️ Additional flap closure tapes which shall be used to close the flaps for further protection of the zip area.
- 🛡️ Made with finest needles of 13mm and the seams are sealed with seam sealing tapes by heat transfer method.

MRP
₹1,850/-





SHOE COVER x 1 PAIR

Maximum Height for better protection.



EYE GOGGLES X1

Complete coverage, end to end protection.



FACE SHIELD X 1

Anti fog treated, Foam for comfortable headrest, Elastic fastener to fit every head size.



NITRILE GLOVES x 1 PAIR

High quality powder free nitrile gloves.



3 PLY MASK X 2

3 Layers with anti bacterial filter and nose pin for better comfort.



BIODEGRADABLE DISPOSABLE BAG X 1

All the items above after use can be collected in this bio hazard bag and discarded.



**COVID19/INMAS/
140520-17/RATL-01**

ZEB-FM500G

Disinfectant Fogging Machine

MRP
₹8,999/-




Air vent for optimized airflow



Fog volume regulator



Power switch

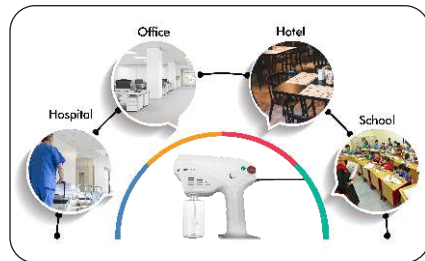

User guide

Step 1



Connect the Power cord of the fogging machine to the power outlet & turn on the red switch present on the left side of the machine to warm-up.

Step 2



Once the switch is on warm-up, the fogging machine starts to Preheat. Once it is heated. The LED light will start to glow. (Pre-heating takes about 40secs)

Step 3



Remove and fill the container with a mixture of sanitizer & water. Please make sure the mixture is 1:1 ratio. (Kindly note that the sanitizer used should not be in the form of gel)

Step 4



Fix the container back into the machine.

Step 5



Rotate the knob (Fog volume controller) to the right as required. The fogging machine will now start to produce fog.

**HANDY
DESIGN**

**EASY
OPERATION**

**200ML
CAPACITY**

**3 METER
CABLE**

After use

Step 1



Rotate the knob (Fog volume controller) to the extreme left.

Step 2



Remove the container from the fogging machine.

Step 3



Again, Rotate the knob (Fog volume controller) to the extreme right. Now the fogging machine will remove all the excess sanitizer from the tube. Rotate back the knob to the extreme left as soon as the machine stops producing any fog.

Specification

Input voltage	220VAC
Power consumption (max)	1200W
Tank capacity	200ml
Housing material	ABS
Power cord length	3 meters
Product Dimension (W x D x H)	220 x 120 x 210 mm
Net. Weight	1.3 kg

Caution

Do not use any liquid with oil

Keep away from children

Fire risk! Do not use any disinfectant with alcohol content

Following alcohol free liquids can be used.

- | | |
|----------------|----------------|
| 1. Glycerol | 5. Nano-silver |
| 2. Glyol | 6. Essence |
| 3. Glycerinum | 7. Water |
| 4. Bactericide | |



Disinfectant Fogging Machine

MRP
₹8,999/-



User guide

Step 1



Connect the Power cord of the fogging machine to the power outlet & turn on the switch present on the rear side of the machine to warm-up.

Step 2



Once the switch is on warm-up, the fogging machine starts to Preheat. Once it is heated. The LED on the top will start to glow. (Pre-heating takes about 2-3 mins)

Step 3



Take a container with a mixture of sanitizer & water. Please make sure the mixture is 1:1 ratio. (Kindly note that the sanitizer used should not be in the form of gel.)

Step 4



The container to be placed in the holder on the fogging machine & the intake tube is to be inserted to the container.

Step 5



Rotate the silver knob (Fog volume controller) to the right as required. The fogging machine will now start to produce fog.

**PORTABLE
DESIGN**

**2000Cu.ft/min
OUTPUT**

**LED
INDICATOR**

After use

Step 1



Rotate the knob (Fog volume controller) to the extreme left.

Step 2



Remove the container from the fogging machine.

Step 3



Again, Rotate the knob (Fog volume controller) to the extreme right. Now the fogging machine will remove all the excess sanitizer from the tube. Rotate back the knob to the extreme left as soon as the machine stops producing any fog.

Features

- Large coverage in less time.
- Easy to carry
- Safe & User friendly

Specification

Input voltage	110Vac-240Vac
Power consumption (max)	1500W
Fogging output	2000Cu.ft/Min
Spray distance	6 meters
Housing material	Plastic
Power cord length	1 meters
Product dimension (W X D x H)	300 x 300 x 200 mm
Package dimension (W X D x H)	385 x 385 x 235 mm
Net.weight	2.5kg

Caution

- Keep away from children
- Fire risk! Do not use any disinfectant with alcohol content
- Following alcohol free liquids can be used.
 1. Glycerol
 2. Glycol
 3. Glycerinum
 4. Bactericide
 5. Nano-silver
 6. Essence
 7. Water

Package contents

Fogging machine	1 Unit
User manual	1 Unit
Power cord	1 Unit

MRP
₹7,999/-

AD801 Contactless Infrared Thermometer



Contactless scanning

Safe and hygienic way to check health status in you homes, offices, malls, factories etc. Easy to use device with large LCD display.

Normal



Fever



High Fever



High Accuracy



Celsius and Fahrenheit Mode



High Temperature Indication



Large Backlight Display



120g approx(w/ battery)
110g approx(w/o battery)



Compact Size
43x77x149mm



Rapid Operation



Auto Shutoff
For Power Saving



2 x AAA Battery



16 Set Memory



ZEB-DT100

Contactless
IR Thermometer



MRP
₹6,499/-

Vocal for Local
#Atmanirbhar Bharat



Contactless scanning

Safe and hygienic way to check health status in your homes, offices, malls, factories etc. Easy to use device with large LED display.



Sleek Lightweight Design



High Accuracy



Celsius & Fahrenheit Mode



High Temperature Indication



Large LED Display



66g approx(w/ battery)
50g approx(w/o battery)



Compact Size
37 x 37 x 138mm



Rapid Operation



Auto Shutoff
For Power Saving



2 x AAA size
Battery



50 Set
Memory

Features

Contactless Infrared Thermometer

Easy to use device with Large LED Display

High Temperature Indication

High Accuracy & Rapid Operation

Specification

Measurement mode Non-Contact

Display size H 26mm, W 21mm

Range of measurement 32°C - 43°C

Measuring distance 1.5cm - 5cm

Measuring time 1s

Accuracy

±0.2°C

Memory data

50 sets

Relative humidity

≤80%

Auto Shutdown

16s

Operating ambient temperature

16°C - 35°C

Power Source

DC 3V (AAA battery x 2)

Product Dimensions (W x D x H)

37 x 37 x 138 mm

Product Weight (without battery)

50g

Package contents

IR Thermometer

1 Unit

User manual

1 Unit

Multifunction UV Sterilizer Box

ZEB-UVS500

99.9% Sterilization
Rate


Dual row UV ray lamps produce 273.7nm band ultraviolet rays, which cuts off DNA/RNA from the micro organisms and kills them.



MRP
₹2,499/-



Effective against
CORONA VIRUS
COVID-19

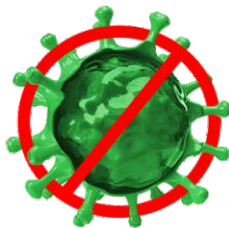

Aromatherapy

Add few drops of any essential/perfumed oil in the slot, then press the incense button.



More than just smartphones

Suitable for a variety of items like masks, handkerchieves, watches, smartphones, keys, glasses, earphones, hearing aids, makeup brushes, nail cutters etc.



Easy to use

Simple and easy to use with notification LEDs.



UV
Disinfection



Intelligent
Sterilization



Quiet
Operation



Spacious
Design



Type C
Input



Dual USB
Output



Aroma
Therapy

Features

99.9% sterilizing rate

360° sterilization in just 3 minutes

Suitable for masks, smartphone, accessories, glasses, etc

Can add Aroma*

Dual row UV Lamps

LED Indicators

Specification

Ultraviolet wavelength	253.7nm
Input power	DC 5V 1A
Power consumption (Diffuse)	1W

Power consumption (Sterilize)	2W
Max power consumption	5W
Supported phone size	upto 16.5cm
No. of USB Output	2 (DC 5V 1A max)
Inner cavity dimension	104 x 180 x 30mm (W x D x H)
Product dimension	125 x 232 x 51mm (W x D x H)
Package dimension	145 x 65 x 260mm (W x D x H)
Net. Weight	306g

Package contents

ZEB-UVS500	1 Unit
Cable	1 Unit
User manual	1 Unit

*Needs few drops of any essential/perfumed oil.

ZEB-UVS200L

Portable UV Sterilizer

99.9%
Sterilizing Rate



MRP
₹1,999/-



ZEB-UVS200L is a portable UV sterilizer that makes your life easier with a simple press and sterilize operation. It eliminates 99.9% of germs, viruses and bacteria on any surface without chemicals or irritants. It can also destroy odour causing bacteria and mold.



ZEB-UVS200L comes in compact and foldable design. It can be easily carried in a hand bag, back pack or luggage to disinfect. It can be used on various objects like keyboard, laptop, mouse, mobile phone, toys, cloths, comb, keys and more.



UVC Light



Portable



LED Indicator



Foldable Design

Features

Kills upto 99.9% of germs on hard surfaces

Sanitize without chemicals or irritant

Sterilization in just ~30 seconds

Destroys odor-causing bacteria

Suitable for keyboard, laptop, mouse, phones, toys, comb etc

Specification

Ultraviolet Wavelength	254nm
Input Power	DC 5V 1.5A (AAA x 4)
Power consumption (Diffuse)	1W
Power consumption (Sterilize)	2W
Max Power consumption	2W
Bulb Life	20,000 H

Operating Temperature	0-50°C
Product Dimension (W x D x H)	125 x 35 x 28mm
Net. Weight	65g

Package contents

ZEB-UVS200L	1 Unit
User Manual	1 Unit

Caution

UV radiation Hazard!
Protect your eyes and skin from exposure.
Keep away from children

Fingertip Pulse Oximeter ZEB-FPO500

MRP
₹4,499/-

CE/FCC/RoHS



ZEB-FPO500 Fingertip Pulse Oximeter is a quick and precise way to check pulse rates and blood oxygen saturation levels. The finger clamp plus simple one-button design allows for easy operation. Small portable size makes it easy to handle and carry. It can be used at home or hospital for better pulse rate monitoring.



Heart Rate
(BPM)



Blood-Oxygen
(SpO2)



Energy
Efficient



One Key
Control



Low-Battery
Indication



Perfusion
Index



The advanced sensors inside the ZEB-FPO500 fingertip pulse oximeter gives you SpO2 level and pulse rate in a very compact form factor.

The ZEB-FPO500 fingertip pulse oximeter is a very easy to use device. It has very simple operation and easy to read digital display.



Features	
Heart rate measurement	
Blood-Oxygen measurement	
Auto shutdown	
Low battery indication	
LCD Display	
Compact design	
High coverage finger pad	
Specification	
Display type	LCD Display
Screen size	14 x 28mm
SpO2	
Measuring range	0-99%
Resolution	1%
Accuracy	70-99%, ±2%; 0-69%, no definition

Pulse rate	
Measuring range	25-250BPM
Resolution	1 BPM
Accuracy	± 1BPM/± 1% BPM (whichever is the larger value)
Input power	2 x 1.5V AAA Batteries
Housing material	ABS
Operating temperature	5°C - 40°C
Product dimension (W X D x H)	32 x 58 x 34 mm
Package dimension (W X D x H)	86 x 70 x 40 mm
Net.weight (w/o batteries)	26.5g

Package contents	
ZEB-FPO500	1 Unit
User manual	1 Unit

ZEB-AM500P

Portable Multi-purpose Atomizer



MRP
₹599/-

ZEB-AM500P is a handy multipurpose atomizer that can be used for everyday care. It comes in small and light weight design, which can be easily carried in the backpack or purse. It can be used to sanitize with 75% surgical alcohol or used as a water mist sprayer to refresh your skin by keeping it hydrated.



Hand sanitizer spray with 75% surgical alcohol (ethanol for disinfection)



Facial mist steamer with water



Rechargeable
Battery



30ml
Capacity



Mini Mist
Sprayer



Ultra
Portable



Stylish and Compact

ZEB-AM500P comes in a stylish and modern design making it more aesthetic. Its compact size makes it to fit anywhere according to your convenience.

Colors

 Pink

 White

 Green



Features

- Can be used as facial steamer with only water
- Moisturize and refresh your skin
- Easy to use
- 30ml capacity
- Compact and handy design
- Rechargeable battery with protection board

Specification

Fogging	1.25 - 1.45ml/60s
Capacity	30ml
Input voltage	DC 5V
Battery	200mAh
Charging time	2 H
Material	ABS + PC
Product dimension (W X D x H)	36 x 100 mm
Package dimension (W X D x H)	41 x 30 x 120 mm
Net.weight	48g

Caution

- Use only water or 75% surgical alcohol (ethanol for disinfection).
- Do not use any liquid with oil.
- Do not use on face if alcohol is added.
- Fire Hazard ! Do not use near fire.
- Keep away from children.

Package contents

ZEB-AM500P	1 Unit
Charging cable	1 Unit

V0620



- India's most admired IT Peripherals & Gaming Brand **WPI 2019**
- Promising Brand **The Economic Times 2019**
- Best Brand for Corporate Gifting (Gadgets) **DT 2018**
- Best Indian IT Brand of the year 2018 **NCN 2018**
- Brand of the year (Speakers) **Silicon India 2018**

All brand names, logos and registered trademarks mentioned are property of their respective owners.

IT Peripherals • Sound Systems • Surveillance
Mobile & Lifestyle Accessories • Power Solutions
Health Care Products

128+
Service Centers

www.zebronics.com

