




Electric Toothbrush

ZEB DENTL20

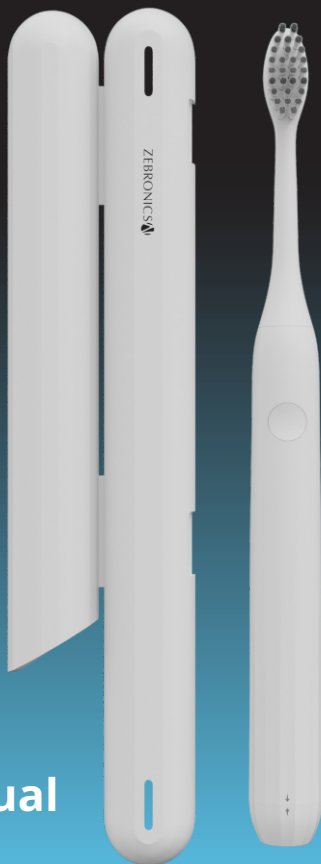

User Manual

www.zebronics.com

Dear Customer,
Please refer to the user manual for optimum usage of **ZEB-DENT L20** Electric toothbrush.

Features:

- Up to 30000 Vibrations per minute - Sonic Technology tooth brush
- Smart Quad Timer
- IPX7 Rating
- 3 Modes (Sensitive, Clean, White)
- Charge lasts up to 40Days*
(2 Mins x 2 Times a day)
- Travel Case included
- 2 x Additional Heads

Specifications:

Input Voltage	: DC 5V
Charging time	: 1 Hours approx
Material	: ABS Plastic
No of Brush heads	: 1 + 2 (Additional Brush heads)
Usage	: 40 Days* (2 Mins x 2 Times a day)
Product dimension (Dia x H)	: 1.85 x 20.5 in cm
Net Weight	: 37g

*Disclaimer : *Backup for standard usage of daily 2 times for 2 Minutes.
Actual backup time might vary depending upon usage."*

Package Content

Electric Toothbrush	: 1 U
Additional Heads	: 2 U
QR Code Guide	: 1 U
Travel Case	: 1 U
Charging Cable	: 1 U

Charging the tooth brush:


- At the bottom of the tooth brush there is a detachable cap - which is indicated with the arrow marks as mentioned in the image below



- Turn the bottom lid counter clockwise - the cap is now unlocked - please pull the cap and it will come off the charging port.



- Now charge the tooth brush with any bundled type-C charging cable from DC 5V rated adapter rated 5V
- Kindly close the charging port as soon as the charging is done, also the charging port cap is one of the vital reasons for the product to be safe from water log.
- For closing the charging port protection cap, align the cap with the tooth brush body in a way that the arrow mark in the inside of the cap matches with the arrow mark of the body.

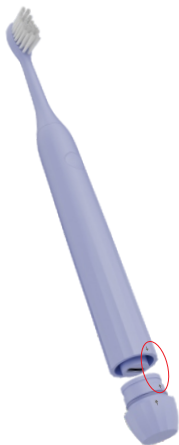

- Now press the cap so that the cap gets inserted into its cavity and looks like this.



- Now turn the bottom cap clockwise till the arrow marks align together.

Usage Instructions:

- This model comes with a single press button for operation. Single press for switching the toothbrush on.

- After Switching on, use a single press for choosing the desired speed mode - Clean, White and Sensitive.
- After choosing the desired speed mode, use this brush like normal toothbrush to clean one quadrant of the mouth
- For every 30 seconds, there will be an indication to switch the Quadrant.
- After the completion of 4 quadrants for 30 seconds per Quadrant, the brush will turn off automatically.
- In case to stop the tooth brush mid way anytime after the first 10 seconds, single press the button to stop the tooth brush.

Cleaning the tooth brush:

Just like a normal tooth brush, this toothbrush can be rinsed with normal free flowing water. Strictly the usage of any solvents, acids, etc to be avoided.

Changing the brushing heads:

- To change the brushing heads, pull the brushing head from the tooth brush body.



- Similarly to connect the bristle head to the body of the tooth brush in a way that the bristle head when facing upward should align with the on / off button.

ISO 9001:2015 Certified Company
www.zebronics.com

